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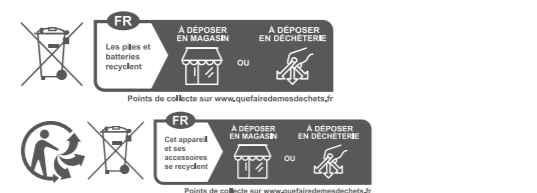
UK REP

EVATOST CONSULTING LTD
Suite 11, First Floor, Moy Road Business Centre, Tall Walk, Carriff, Wexley, CP15 10N
contact@evatmaster.com

EC REP

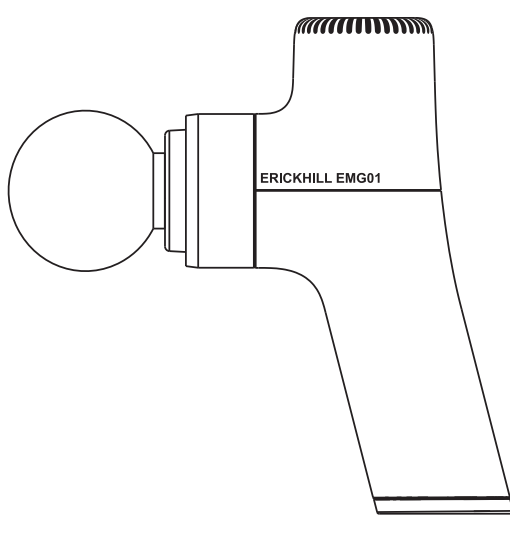
Evatmaster Consulting GmbH
Bismarckstr. 30
60323 Frankfurt am Main, Germany
contact@evatmaster.com

Erickhill

User's Manual

Massage Gun **EMG01**

CE RoHS FC ERII
MADE IN CHINA
Contact us: support@erick-hill.com

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UK REP

EVATOST CONSULTING LTD
Suite 11, First Floor, Moy Road Business Centre, Tall Walk, Carriff, Wexley, CP15 10N
contact@evatmaster.com

EC REP

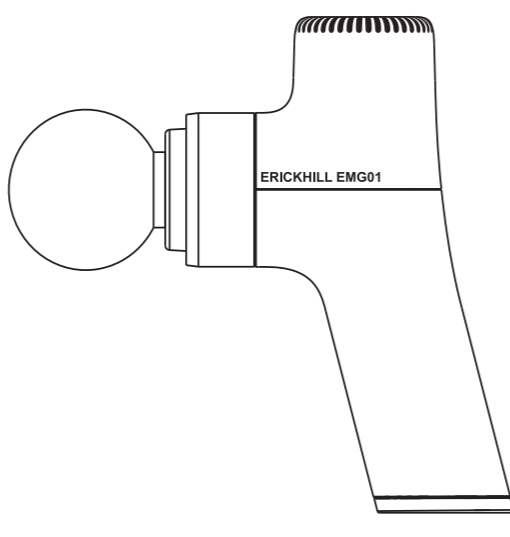
Evatmaster Consulting GmbH
Bismarckstr. 30
60323 Frankfurt am Main, Germany
contact@evatmaster.com

Erickhill

User's Manual

Massage Gun **EMG01**

CE RoHS FC ERII
MADE IN CHINA
Contact us: support@erick-hill.com

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Notes

Warnings

- Do not use this product on sensitive areas such as neck and head.
- Do not use this product directly on thoracic, lumbar spine, and other areas of nerve concentration.
- Do not use this product on areas such as the face, back of hands, etc.
- Do not use this product on protruding bones and joints.
- Do not use this product when you have obvious skin, muscle, or joint injuries.
- Do not use this product when there are already signs of pain or other discomfort.
- Children under 8 years old are prohibited from using this product. Children are not allowed to use, clean, or maintain this product without supervision.
- Do not place any objects in this product's air vents.
- Do not immerse this product in water, and do not allow water to enter this product's air vents.
- Do not put the product into fire or near fire, high temperatures and direct sunlight.
- Do not use other charging cables, please use the special charging cable corresponding to the model of this machine when charging, and each charging should not exceed 3 hours.
- Please do not disassemble or modify it in any way.
- Do not touch the joint between the body and the massage head while the device is operating (with or without the massage head installed) to prevent pinching.
- Do not charge for a long time, please unplug the charging cable when the power prompt is full.
- Please dispose of this device, battery and other accessories according to local regulations and do not dispose of them as household waste. Improper handling of the battery may cause the battery to explode.
- Do not touch when the device starts.

Recycling

When the service life of this product is over, please do not throw it away together with general household waste, put this product in the government-designated recycling station, or send this product back to the manufacturer's after-sales service center, and the personnel of the service center will use environmentally friendly way of processing.

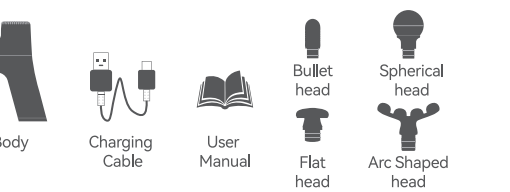
Statement

If any of the following apply to you, please consult a doctor before use and follow the doctor's guidance in using this product.

- Pregnant.
- Using electronic medical devices such as pacemakers.
- Recently had surgery.

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Product introduction

Body Charging Cable User Manual Flat head Arc-Shaped head

Usage

Power Button/Frequency gear switch Button
Power On/Off
Long Press for 2 seconds to turn on or off the device in standby mode. It will automatically shut down if there is no operation for 15 seconds.
Frequency gear switching
In the power-on state, short press the button to switch the frequency gear.
1st-Gear-2000r/min, 2nd-2600r/min, 3rd-3200r/min.
It will automatically enter the standby state after 10 minutes of continuous use, (automatically shut down if there is no operation for 15 seconds), and it needs to be restarted when it is used again.
After the massage head is blocked, it stops outputting, and it is set for the protection function of abnormal operation, from running to standby mode, short press the button to restart.

Battery-Charging Indicator
Blue light on: battery remaining 51%-100%.
Purple light on: battery remaining 25%-4%.
Red light flashing: less than 1%, and it needs to be charged in time.

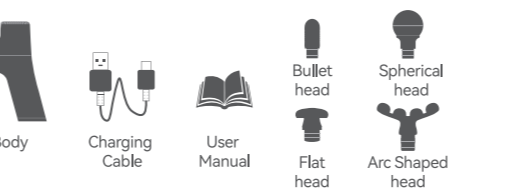
USB type-C charging port
Before using the product for the first time, please fully charge it.

When charging, charging indicator
Red light breathing: battery 1%-4%.
Purple light breathing: battery 5%-50%.
Blue light breathing: battery 51%-95%.
Blue light on: battery 96%-100%.

CE RoHS FC ERII
MADE IN CHINA
Contact us: support@erick-hill.com

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Product introduction

Body Charging Cable User Manual Flat head Arc-Shaped head

Usage

Power Button/Frequency gear switch Button
Power On/Off
Long Press for 2 seconds to turn on or off the device in standby mode. It will automatically shut down if there is no operation for 15 seconds.
Frequency gear switching
In the power-on state, short press the button to switch the frequency gear.
1st-Gear-2000r/min, 2nd-2600r/min, 3rd-3200r/min.
It will automatically enter the standby state after 10 minutes of continuous use, (automatically shut down if there is no operation for 15 seconds), and it needs to be restarted when it is used again.
After the massage head is blocked, it stops outputting, and it is set for the protection function of abnormal operation, from running to standby mode, short press the button to restart.

Battery-Charging Indicator
Blue light on: battery remaining 51%-100%.
Purple light on: battery remaining 25%-4%.
Red light flashing: less than 1%, and it needs to be charged in time.

USB type-C charging port
Before using the product for the first time, please fully charge it.

When charging, charging indicator
Red light breathing: battery 1%-4%.
Purple light breathing: battery 5%-50%.
Blue light breathing: battery 51%-95%.
Blue light on: battery 96%-100%.

CE RoHS FC ERII
MADE IN CHINA
Contact us: support@erick-hill.com

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Product introduction

Body Charging Cable User Manual Flat head Arc-Shaped head

Usage

Power Button/Frequency gear switch Button
Power On/Off
Long Press for 2 seconds to turn on or off the device in standby mode. It will automatically shut down if there is no operation for 15 seconds.
Frequency gear switching
In the power-on state, short press the button to switch the frequency gear.
1st-Gear-2000r/min, 2nd-2600r/min, 3rd-3200r/min.
It will automatically enter the standby state after 10 minutes of continuous use, (automatically shut down if there is no operation for 15 seconds), and it needs to be restarted when it is used again.
After the massage head is blocked, it stops outputting, and it is set for the protection function of abnormal operation, from running to standby mode, short press the button to restart.

Battery-Charging Indicator
Blue light on: battery remaining 51%-100%.
Purple light on: battery remaining 25%-4%.
Red light flashing: less than 1%, and it needs to be charged in time.

USB type-C charging port
Before using the product for the first time, please fully charge it.

When charging, charging indicator
Red light breathing: battery 1%-4%.
Purple light breathing: battery 5%-50%.
Blue light breathing: battery 51%-95%.
Blue light on: battery 96%-100%.

CE RoHS FC ERII
MADE IN CHINA
Contact us: support@erick-hill.com

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Product introduction

Body Charging Cable User Manual Flat head Arc-Shaped head

Usage

Power Button/Frequency gear switch Button
Power On/Off
Long Press for 2 seconds to turn on or off the device in standby mode. It will automatically shut down if there is no operation for 15 seconds.
Frequency gear switching
In the power-on state, short press the button to switch the frequency gear.
1st-Gear-2000r/min, 2nd-2600r/min, 3rd-3200r/min.
It will automatically enter the standby state after 10 minutes of continuous use, (automatically shut down if there is no operation for 15 seconds), and it needs to be restarted when it is used again.
After the massage head is blocked, it stops outputting, and it is set for the protection function of abnormal operation, from running to standby mode, short press the button to restart.

Battery-Charging Indicator
Blue light on: battery remaining 51%-100%.
Purple light on: battery remaining 25%-4%.
Red light flashing: less than 1%, and it needs to be charged in time.

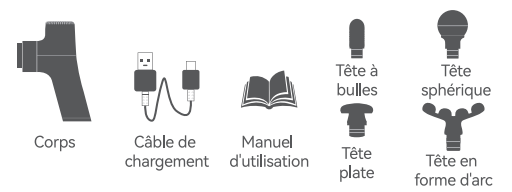
USB type-C charging port
Before using the product for the first time, please fully charge it.

When charging, charging indicator
Red light breathing: battery 1%-4%.
Purple light breathing: battery 5%-50%.
Blue light breathing: battery 51%-95%.
Blue light on: battery 96%-100%.

CE RoHS FC ERII
MADE IN CHINA
Contact us: support@erick-hill.com

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